



MERIVAARA TROLLEY

Antistatic mattresses

Title: Merivaara Antistatic Trolley Mattress
Codes 51751 / 51781

Purpose of use


- Mattresses for trolleys in hospitals.

Product structure

- Core R50200 (50 kg/m³) fire retardant HR-foam.
- Cover black antistatic polyurethane coated hygienic fabric, which is impermeable to liquids but still breathable. Zipper on the end.
- Velcros on the bottom and 4 pcs on the sides for fixing. 51781 mattress has also 2 bands on the sides.

Benefits

- Secure and easy-care mattresses.
- With several velcros mattress is easy to fasten on the trolley.

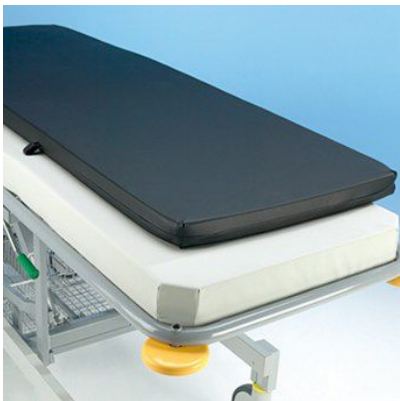

Cleaning and maintenance

- Disinfection by wiping according to the instructions with diluted chlorine-amine or an alcohol-based detergent or ethanol /chlorine max 1000 ppm, rinse and dry thoroughly. If necessary, wash the cover 70-95°C inside out, gentle tumble drying at low temperature.
- Replace broken cover to a new one.

Flame retardancy & CE

Flammability by standard IMO FTPC part 9. Merivaara Trolley mattresses are CE-marked. Country of origin Finland.

Size

Standard sizes 65*193*5 and 65*193*8 cm